


T A O clean


Sonic Toothbrush

& UV Sanitizing Base

Instruction Manual

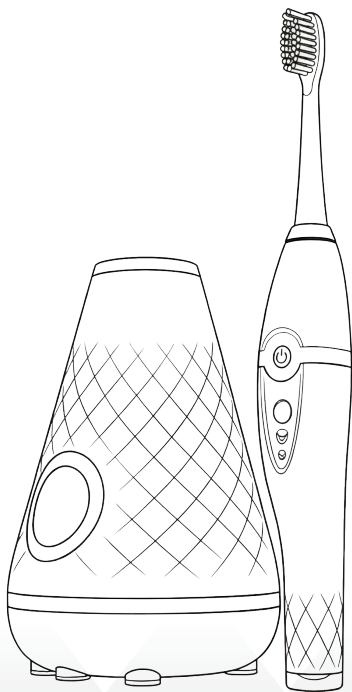


WHAT'S IN THE BOX

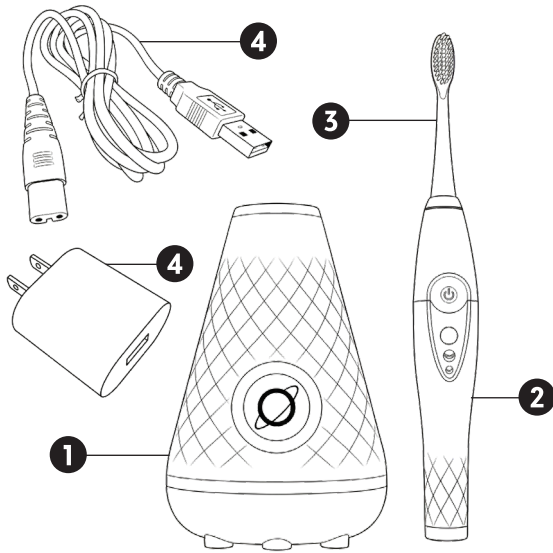

- ❶ UV Sanitizing Base
- ❷ Toothbrush Handle
- ❸ Daily Care Whitening Toothbrush Head
- ❹ Charging Cord and Power Brick

Contents may be hidden under cardboard within the box. Open the flap to find all components.

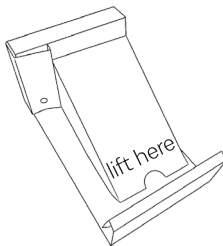


GETTING STARTED

1. Setup the UV Base

Plug In: Connect the charging cord to the back of the UV Sanitizing Base and the power brick.

Power On: Plug into a wall outlet; the base should light up.

Stay Connected: Keep the base plugged in at all times to ensure the drying, charging, and sanitizing functions work.

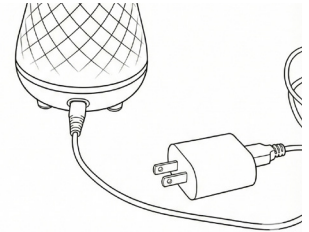

2. Assemble the Brush


Attach Head: Press the brush head firmly onto the handle until secure.

No Twisting: Do not twist the head during attachment or removal.



3. Dock the Brush


Dock the toothbrush upside down in the base. Make sure the toothbrush is docked properly to ensure the drying, charging, and sanitizing functions are activated.

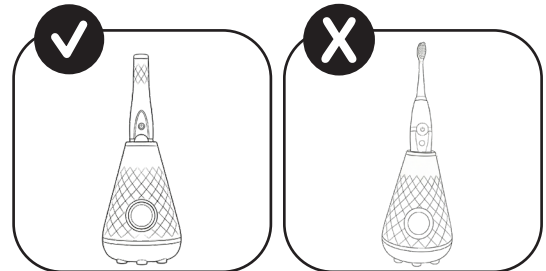

BRUSHING BASICS


Prep: Apply toothpaste and place the brush in your mouth before powering on to prevent splashing.

Power & Modes: Press the power button to start. Use the mode button to toggle between **Super Clean** (strong vibration) and **Sensitive Clean** (gentle experience).

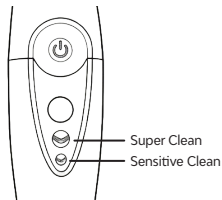

Technique: Gently glide the bristles across teeth and gums; do not use heavy pressure.

Timer: The brush pauses every 30 seconds to prompt moving to a new quadrant of your mouth and shuts off automatically after 2 minutes.

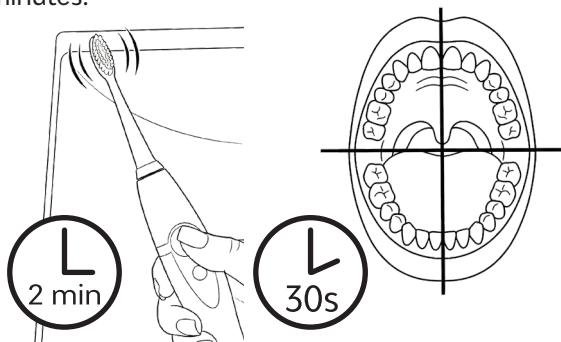


CLEANING & MAINTENANCE


Docking: After brushing, rinse the brush head and shake off excess water.

Proper Fit: Place the handle upside down in the base with the power button facing the logo. This triggers automatic charging, drying, and UV sanitizing.

Base Care: Unplug and wipe the base with a soft cloth once a week to remove residue.

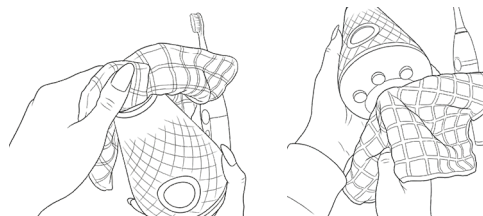

Replace Heads: Change your brush head every 3 months or when bristle color fades. To remove brush heads, grip firmly and pull straight off. Replacements are available on taoclean.com.

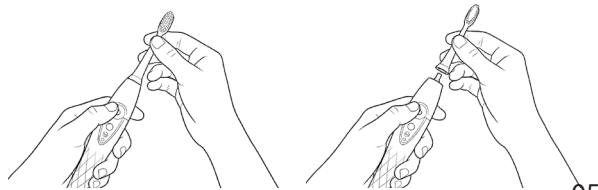

Where can I get replacement brush heads?

Replacement toothbrush heads can be purchased online at taoclean.com/refills. All brush heads ship free, and you can save up to 20% with a brush head subscription.

How do I charge the device correctly?

Charging is automatic. When your toothbrush is docked correctly in the base (see page 3), it triggers a charging, sanitizing, and drying cycle. The handle is made to only fit flush one way and will sit naturally in the correct position.

Why isn't my toothbrush charging?

Try this quick checklist:

- Confirm the charging cord is plugged firmly into the power brick & power brick is plugged into the wall outlet.
- Confirm the charging cord is fully inserted into the base
- Re-dock with your toothbrush head correctly (see page 3)
- Confirm the badge light turns on & pulses blue when your toothbrush is inserted

If it still won't charge, contact support.

Why isn't the logo glowing?

This is normal behavior based on the UV cycle and Night Light settings.

- UV runs for 12 minutes after docking
- The Tao Clean logo pulses blue for the first 5 minutes only
- After that, the badge is either solid white or dark depending on Night Light Mode
- UV runs again 12 hours later, then every 24 hours if not used (timing resets with use)

To change Night Light Mode, press the black button on the bottom of the base.

How long does the battery last?

- Docking daily keeps the battery topped off automatically
- When traveling without the base: the handle will work for 2 to 3 weeks on a full charge (varies by usage)

Note: The base must be plugged into an outlet at all times in order to charge your toothbrush.

CONTACT

Need Help?

For assistance:

Email: support@taoclean.com

Phone: +1-714-363-3599

Web: taoclean.com/pages/start

WARRANTY

TAO Clean toothbrush handles and UV sanitizing bases purchased directly from TAO Clean or from an Authorized Seller are covered by a Limited Lifetime Warranty. Full warranty details are available at taoclean.com

WARNING

UV Exposure: This device contains a UV-C emitter. Do not stare into the light source, as it can be harmful to eyes and skin. Use the device only on flat, horizontal surfaces.

Water & Electricity: If the base falls into water, unplug it immediately.

Care: Do not place any components in the dishwasher.

If the brush head, handle, charger or base become damaged, stop using it.

If you experience pain or discomfort while using, stop and contact your dentist

This equipment has been tested and complies with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Thank you for choosing

T A O clean

Happy Brushing!

Loving your TAO Clean Sonic Toothbrush?
Leave Us a Review!

taoclean.com